CANTEC

- helping people

Fall Sensor

Bruksanvisning | User Guide

Cantec AS

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Read first: Security information

Battery replacement shall be done by authorized persons. The usage is only allowed with recommended battery. CAREFULNESS – danger of explosions may occur if wrong types are used. Used batteries must be disposed according to instruction.

Always read and follow the security information when this symbol is shown.

Important information

In all system that are using radio and network communication there is a potential risk of interference that is out of control for the user.

CareTech's products are constructed to minimize those risks. Despite this you have to be aware that parts of the system can be exposed to different types of disturbance with impaired or intermittent function as a result.

Therefore it's always important to regularly verify all parts in the system, especially radio communication that must be working in all areas. Immediately contact your supplier in suspicions of failure.

The user should be aware of the risk of disturbance from products that communicates within the same or adjacent frequencies.

The Fall Sensor is a highly reliable product but it can't alert in situations that are similar to daily activities and there's also daily activities that are similar to situations when the sensor will alert.

For more information, please contact your supplier or visit www.caretech.se


Package contents

- Fall sensor
- Manual

Description of Fall Sensor

Belt clip

Used for attaching the sensor to the belt.

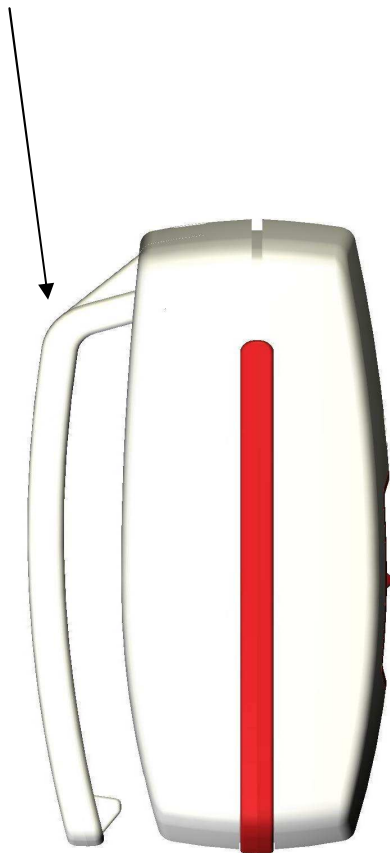Indicator lamp

Red light indicates radio transmission because of fall detection or activated alarm button. Green light indicates acknowledge.

Alarm button

The button is used for manually calling for help.

Functional description

The Fall Sensor contains advanced electronics that register and handles movements in all directions and speed.

In case of rapidly movements the direction and speed are calculated by the sensor to decide if the movements are caused by falling or normal daily activities.

The Fall Sensor also has an alarm button to manually raise an alert. The device is also equipped with automatic battery monitoring and self-testing alarms. If the self-testing alarm fails a silent alert will be sent to the alarm receiver (requires that the function is activated in the carephone).

Installation

The Fall Sensor is compatible with these carephones Gloria+, CareIP and CareIP Mobile. Follow all steps below to activate and install the sensor.

1. Activate Fall Sensor

The Fall Sensor is deactivated by default when delivered and need to be activated before being used. By pressing the red alarm button for approximately 5 seconds the indicator lamp will light green and the unit will be activated. If the sensor is already activated there's no need doing this step.

If you are uncertain or only want to verify the status press the red alarm button for approximately 1 second. Red blink means that the sensor is already activated. No blink means that the sensor needs to be activated.

2. Radio range control

A range control shall always be done at the installation. This will be done by pressing the green button on the carephone until a beep is heard (approximately 4 second). Take a walk in the housing while pressing the alarm button. If the indicator lamp on the Fall Sensor lights red and changes to green then the sensor is in range. If the sensor blinks red it's out of range and the carephone need to be moved (if possible) to improve the range.

3. Programming the Fall Sensor

The Fall Sensor needs to be programmed against a carephone to work. This will be done by pressing the green button on the carephone until a beep is heard (approximately 4 second), push the alarm button on the Fall Sensor and then confirm by pressing the red button on the carephone. For fulfilling the process a test alarm need to be done to verify that the Fall Sensor are communicating against the carephone.

Operation

For the best and most secure function always wear the Fall Sensor in vertical position at waist level as the sensor in this position follows your upper part of the body. Most appropriate attachment is straight ahead by the clip in an ordinary belt or by using our specially designed waist belt. (Not included).

The waist belt is clamped around the waist with the arrow pointing upward. The Fall Sensor is placed with the red alarm button outwards into the small pocket by help of the clip.

The waist belt is an optional part and can be ordered in following sizes: S (95cm), M (100cm), L (110cm), XL (125cm) och XXL (140cm).

When removing the Fall Sensor it's important to handle the sensor with care to avoid unintentional alarm.

Maintenance

Cleaning

The Fall Sensor shall be cleaned by a dampened cloth. No solvents are allowed when cleaning.

The waist belt can be washed in up to 95°C, but in normal case 60°C is good enough. The Velcro shall be put together before washing in order to keep its strength.

Cotton 85% / Lycra 15%

The Fall Sensor must not be exposed for moisture and shall always be removed from the waist belt before washing.

Important! When the Fall Sensor is removed for washing the waist belt it's important that the sensor is placed in correct position when re-inserting the sensor in the waist belt, see operation.

Battery replacement

A signal for low battery is transmitted automatically to the carephone when the battery should be changed. The battery should always be replaced by an authorized person.

See section "Read first: Security information" on page 4 for more information.

Technical data

Dimensions:	52 x 33 x 24mm (HxWxD)
Battery:	3V lithium battery CR2450
Battery life:	Up to 1 year
Communication:	Radio, 869,2125MHz
Water resistance:	According to IP65 (Water jet)

Equipment class

Class 1 radio equipment

Environment

This product is intended for indoor use in a normal residential environment

Temperature:	Operating temperature +5 to +35°C
Humidity:	0 till 75 % relative humidity (non-condensed)
Environmental class:	1

Environmental information

This product complies with the requirements of the EU directives, 2002/96/EC (WEEE).

WEEE regulate the product liability for electrical and battery recycling with the purpose of increasing recycling and minimising waste. The unit is marked with the “crossed out wheeled bin” logo, which indicates that it shall be handed in for recycling.

The product can be returned free of charge to a recycling station that is connected, directly or via a recycling system, to Care Tech or to your distributor. For detailed instructions please check with your distributor or visit our website, www.caretech.se

NOTE! This information and recycling instructions applies to European Union member states only. For other countries please check local legislation or contact your distributor.

Manufactured in accordance with the EU directive, 2002/95/EC (RoHS).

The materials used in the waist belt meet the textile safety of Oeko-Tex standard 100.

The information contained here is subject to change without notice. The only warranties for CareTech products and services are set forth in the express warranty statements accompanying such products and service. Nothing herein should be constructed as constituting an additional warranty. CareTech shall not be liable for technical or editorial errors or omissions contained herein.

www.caretech.se