
## Item list

- 1) Manual
- 2) Control Unit
- 3) Absorgel Pouch
- 4) Power supply (main inlet cables)
- 5) USB Memory Stick
- 6) Time Card
- 7) Padding
- 8) Seal
- 9) Positioner
- 10) Strap
- 11) Pressure Chamber

## Boot Setup first time use

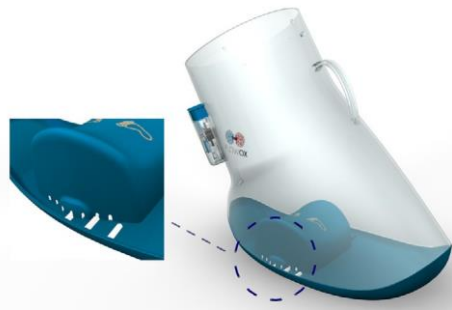1) Place positioner grooves onto Pressure Chamber ridges. Move forwards or backwards as required.

2) Insert Padding, place Air Vents into holes, attach fastening band at the top to the studs.

3) Place Seal onto stud and around the top of the Pressure Chamber.

4) Attach the strap around the bottom of the seal at the front and fasten onto the stud.

## Treatment - Regular use

- 1) Connect the Hose from the Pressure Chamber to the Control Unit.
- 2) Plug in external power supply to the Control Unit before plugging it into the wall socket.

- 3) Gently pull the seal over the outside of the Pressure Chamber.


- 4) Place foot on top of positioner centrally in Pressure Chamber.

- 5) Pull the seal onto the leg to make the Pressure Chamber air tight.

- 6) Start the treatment by switching on the device by pressing the on/off button.

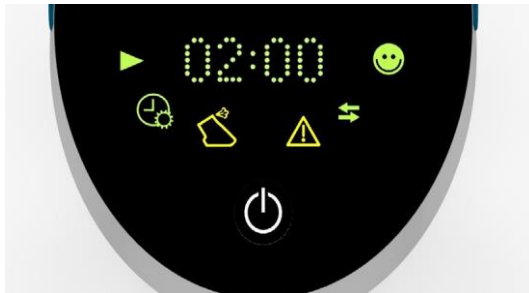- 7) After the treatment turn off the device by pressing the on/off button.

- 8) Roll down the seal onto the outside of the Pressure Chamber.

- 9) Remove the air from the padding by pressing the release valve button.

- 10) To remove air quickly press down on inner Padding.

**E2 Major Air Leakage**  
Higher leakage than normal and treatment not being delivered  
Solution: Check the Seal, Hose with Filter and adjust if necessary. Restart the device by pressing On/Off Button.

**E3 Low Air Volume**  
Not enough air to function (Low air volume)  
Solution: check for obstruction, e.g. kinked hose.  
Restart the device by pressing On/Off Button.

**E4 Data transfer error**  
Data not transferred, e.g. due to taking out USB Memory Stick too early or software error  
Solution: Switch off the device and try inserting the USB Memory Stick into the Control Unit again.

**Time Card** | When available time on goes down to 25 hour or less, hourglass and the time left will show on display  
Solution: Contact your clinician to load more time with Time Card.

**Note:**  
If the problem persists, contact clinician or local distributor.  
The Control Unit can be reset by pressing On/Off button for 13 seconds.  
More information is available in the manual.

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